**oozz**  
**OZ-3611**


Importer: Elite Capital Limited

Address: Studio 4 Clarks Courtyard, 145 Granville Street,  
Birmingham, B1 1SB, United Kingdom.

Website: [oozz.com](http://oozz.com)

Email: [care@oozz.com](mailto:care@oozz.com)

OZ-3611

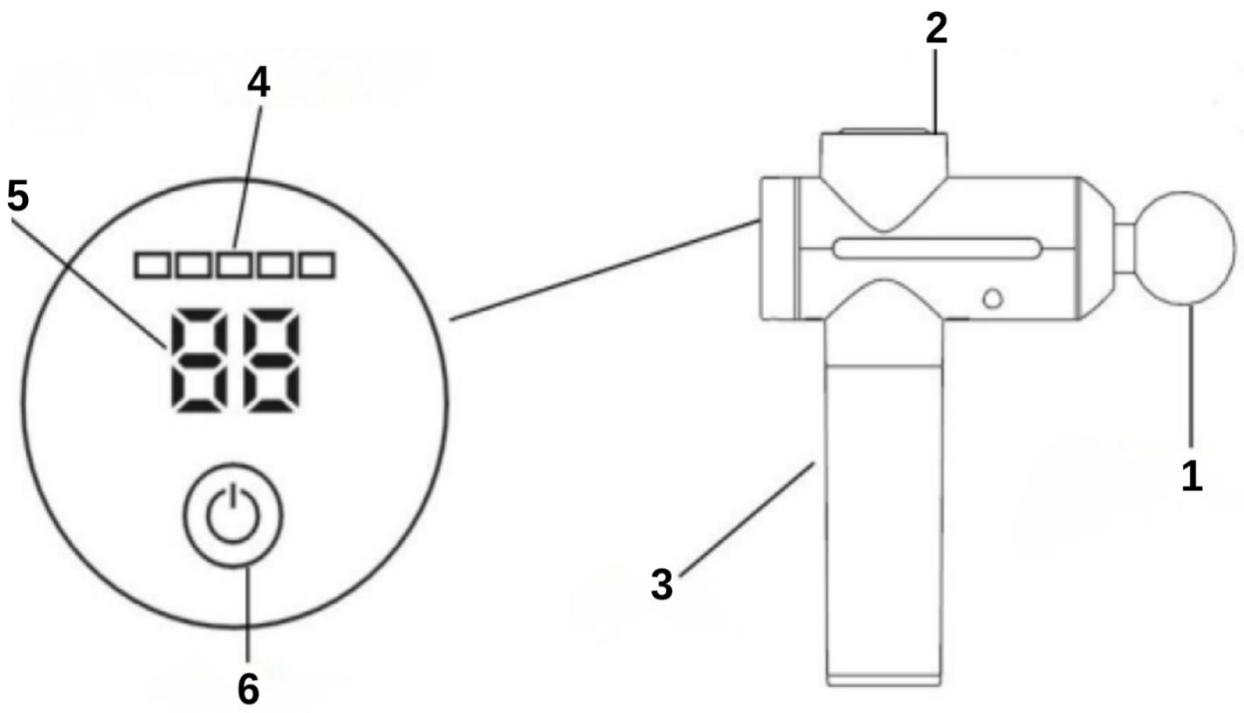7

8

9

10

11

12

13

14

15

Welcome to OOZZ and congratulations on your new purchase! Secure your peace of mind by activating your 12-month warranty and claiming extra 6 months for free at [www.oozz.com/warranty](http://www.oozz.com/warranty). Beyond our warranty, you are also covered by your country's legal consumer rights, ensuring complete protection for your purchase.

1. Massage Head
2. Force Indicator Light
3. Handle
4. Battery Indicator
5. Speed Level Display
6. Power Button
7. Charging Cable
8. Flat Head (Suitable for large muscle groups such as the back, waist, arms)
9. Sphere Head (Ideal for large muscle groups across the body)
10. Cylinder Head (For joints and targeted deep-tissue areas)
11. U-Head (For the spine and area behind the ankle)
12. Silicone Swamm Head (Ideal for the whole body, especially medium to large muscle areas, due to its soft and gentle surface.)
13. Round Scoop Head (For local and small muscle areas)
14. Camber Head (For the legs, arms, and neck)
15. Warm Head (Massage body with hot compress)

#### **Product Details**

Product Name	OZ-3611 Massage Gun
Model Number	OZ-3611
Charging Input	DC 5V 1A
Charging Time	≤4.5 Hours
Usage Time	120 mins (MAX Level)

#### **Charging Instructions:**

- Make sure the device is switched OFF before charging.
- Charge the product fully before you use it for the first time and after a long period of disuse. For optimal performance, it is best to use the device when fully charged.
- When the device power is low, the Battery Indicator Light will flash to remind you to charge the device.
- When the device is charging, the Battery Indicator Light will blink and will stay solid once fully charged.

#### **How to use:**

- Long press the power button for 3 seconds to Turn ON or Turn OFF the device. Press again to start the device at Level 1 by default.

- Press the power button repeatedly to cycle through 12 intensity levels. The LED display will show the current level.
  - Levels 1-10: Constant speed modes (gentle to strong)
  - Levels 11-12: Variable speed mode: intensity gradually increases and decreases for a wave-like massage effect
- The device features a built-in force sensor:
  - Green light: Optimal pressure
  - Blue light: High pressure - ease up slightly
  - Red light: Too much pressure - reduce immediately for safe use

**Precautions:**

- Charge only with DC 5V 1A certified adapter.
- Do not overcharge or charge near water. Do not expose to water or moisture. This device is not waterproof.
- Do not use if damaged or malfunctioning.
- The product should be used by adults only and keep out of reach of children.
- For your safety, do not use the massage gun while it is charging.
- Use the device only on soft tissues of the body without causing pain or discomfort. Do not use it on head, bones of sensitive parts of the body.
- Stop using the device once there is any sign of pain or discomfort.
- Only the attachments on the recommended muscle groups.
- To protect the bones and muscles, avoid long-term fixed-point massage.
- Do not immerse the device in water or allow water to enter the device.
- Do not use on injured or swollen areas, or if you have open wounds, bruises, or skin conditions.
- Consult your doctor before use if you have any medical conditions, such as heart problems, diabetes, or if you are pregnant.
- Do not use for more than 60 seconds on a single area. Move the device slowly across muscles.
- Do not operate while sleeping, drowsy, or under the influence of drugs or alcohol.
- Use only the provided charger and accessories.
- Stop use immediately if you feel pain, discomfort, dizziness, or notice any abnormal reaction.
- Do not disassemble or attempt to repair the device yourself. Contact customer service for assistance.

**Cleaning:**

- Power off and disconnect before cleaning.
- To clean, wipe the surface of the device with a damp towel and dry it with a soft towel.

**Customer Service**

For assistance, please visit our website: [OOZZ.com](http://OOZZ.com). You can also call or WhatsApp us on 0800 193 6789 (UK), +44800 193 6789 (International) or email us at [care@oozz.com](mailto:care@oozz.com).

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