oozz

OZ-3611


Importer: Elite Capital Limited

Address: Studio 4 Clarks Courtyard, 145 Granville Street,
Birmingham, B1 1SB, United Kingdom.

Website: oozz.com

Email: care@oozz.com

OZ-3611

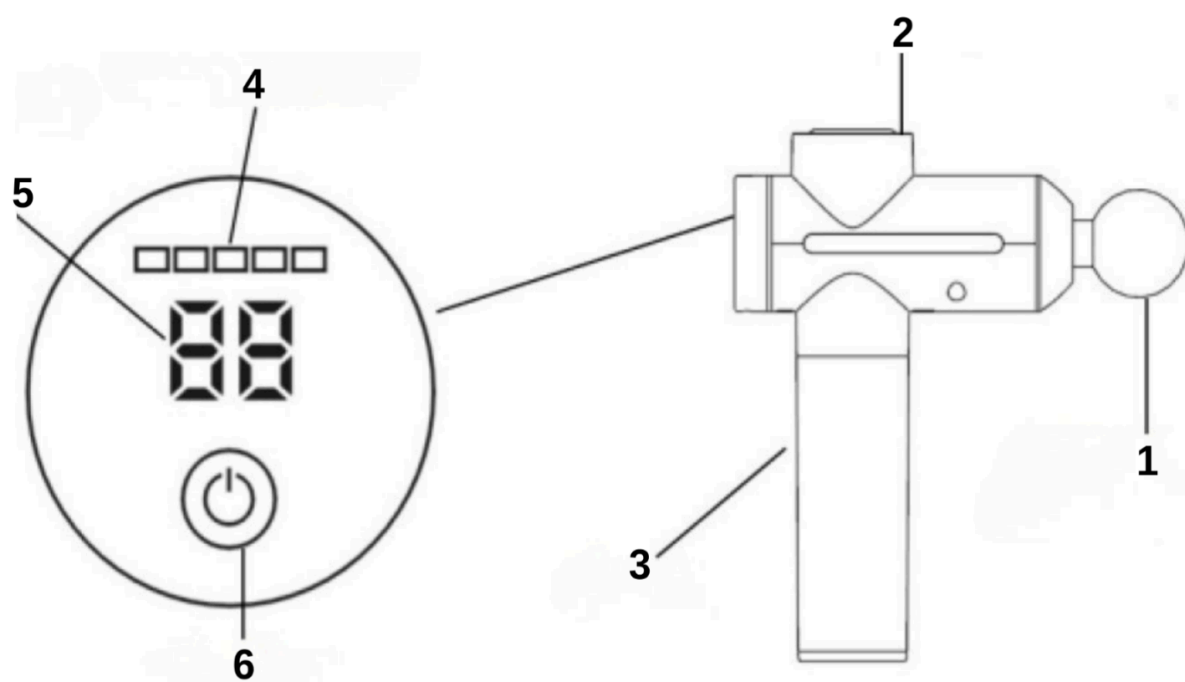7

8

9

10

11

12

13

14

15

Welcome to OOOZ and congratulations on your new purchase! Secure your peace of mind by activating your 12-month warranty and claiming extra 6 months for free at www.oozz.com/warranty. Beyond our warranty, you are also covered by your country's legal consumer rights, ensuring complete protection for your purchase.

1. Massage Head
2. Force Indicator Light
3. Handle
4. Battery Indicator
5. Speed Level Display
6. Power Button
7. Charging Cable
8. Flat Head (Suitable for large flat muscle groups)
9. Sphere Head (Ideal for medium to large muscle groups)
10. Cylinder Head (For deep-tissue massage)
11. U-Head (For around the spine and Achilles tendon)
12. Silicone Swamm Head (For gentler massage of medium to large muscle groups)
13. Round Scoop Head (For small muscle areas)
14. Camber Head (For the legs, arms, and neck)
15. Warm Head (For massaging large muscle areas with a hot compress)

Product Details

Product Name	OZ-3611 Massage Gun
Model Number	OZ-3611
Charging Input	DC 5V 1A
Charging Time	≤4.5 Hours
Usage Time	120 mins (MAX Level)

Massage Heads: Function & Usage

Massage Head Type	Usage
Flat Head	Suitable for large, flat muscle groups in the back, waist or thighs, and aids myofascial release.
Sphere Head	Versatile for use on both large and medium muscle groups.
Cylinder Head	Good for deep-tissue massages on soft areas such as the soles of the feet and palms of the hands.

U-Head	Ideal for massaging around bony areas such as the spine, up to neck, or the Achilles tendon.
Silicone Swamm Head	Suitable for gentler massage of the whole body, especially medium to large muscle areas, due to its soft, springy surface.
Round Scoop Head	For small muscle areas like arms, shoulders, and calves.
Camber Head	Suitable for arms, legs and the neck contours thanks to wider, rounded shape.
Warm Head	Ideal for medium to large muscle areas. Hot compress increases circulation and aids muscle recovery quicker than with standard Massage Heads.

Instructions for Use

1. ON/OFF Button:

- Long press the Power Button to Turn ON the device.
- Press the Power Button once to activate the first speed level.
- Hold the device by the Handle and place the massage head onto the desired body area.
- Long press the Power Button to Turn OFF the device.

2. Inserting and Removing Massage Heads:

- Insert the desired Massage Head into the front opening of the main unit, pressing firmly to the end (Figure 1).
- To remove a Massage Head, gently pull it straight out.

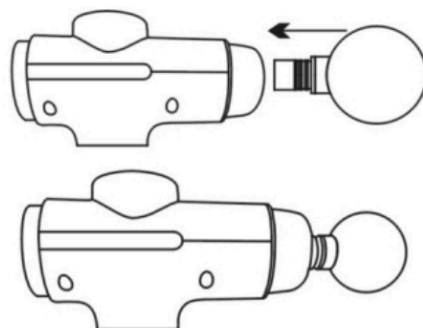
Figure 1

Note: Do not attempt to change Massage Heads while the device is turned on.

3. Speed Settings:

- The Power Button has 12 intensity levels, ranging from 1600RPM to 3200RPM. It starts at level 1, and each press increases the level by one.
- The selected level will be indicated by the number on the LED Speed Level Display.

Note: Regardless of which Speed Level the device is running at when turned off, it will always default back to Level 1 when turned on again.

4. Pressure Indicators:

- The device features a built-in pressure sensor. The current pressure level will be indicated by the Force Indicator Light.
 - Green light: Optimal pressure
 - Blue light: High pressure - ease up slightly
 - Red light: Too much pressure - the device will turn off as a safety precaution

5. How to Charge:

- Connect the provided USB cable directly to the device. Then plug the USB into a 5V 1A USB adapter. The product must be charged via a laptop/computer/power bank USB port. After charging remove the power plug from the USB outlet and then remove the USB out of the product.
- Charge the product fully before you use it for the first time and after a long period of disuse.
- Make sure the device is switched OFF before charging.
- The LED lights on the Charging Interface will light up, confirming the charging process has begun.
- The device's charge level will be shown by however many Indicator lights are solid blue (not flashing):
 - 20% charged: One light
 - 40% charged: Two lights
 - 60% charged: Three lights
 - 80% charged: Four lights
 - 100% charged: Five lights
- A full charge takes about 4.5 hours providing up to 120 minutes of operation. For optimal performance, it's best to use the device when fully charged.

Cleaning & Maintenance

Please note: This appliance is not waterproof so do not immerse the device in water or any other liquid to maintain optimal performance and longevity.

- Before cleaning, ensure the device is switched OFF and disconnected from any charging cables.

- Always use only water or neutral cleaning agents with a damp towel to clean the appliance. Never use caustic or abrasive agents (e.g. vinegar or lime scale removers) since they may cause damage to the appliance. After, dry with a soft towel.
- Keep the appliance away from high temperature and avoid direct sunlight.
- Do not submerge the device or massage heads in water or any other liquid.

Precautions

Charging

- Ensure the device is powered off before connecting the charging cable. Avoid forcefully pulling on the power supply wire or the device.
- Always attach the power cord to the product first, then to the power outlet.
- Charge only with DC 5V 1A certified adapter.
- Do not overcharge or charge near water. Do not expose to water or moisture. This device is not waterproof.
- Do not plug this device into a voltage system that is different from the voltage system specific on the device or charger.
- Prevent overcharging by unplugging the device once the designated charging time is reached.
- Do not charge on soft surfaces like beds and bedding, especially in enclosed spaces. Should the product overheat during charging, disconnect it immediately.
- Avoid charging in high-temperature areas or near open flames to mitigate fire risks.
- Avoid damaging, breaking, or modifying the USB cable. Do not stretch, bend or twist the USB cable excessively, and do not place heavy objects or trap it between objects.
- If the power cord becomes damaged, replace it using only compatible cords or components purchased directly from the original manufacturer.

Usage

- Use only the provided charger and accessories.
- To protect the bones and muscles, avoid long-term fixed-point massage.
- Use the device only on soft tissues of the body without causing pain or discomfort. Do not use it on head, bones of sensitive parts of the body.
- Do not use if damaged or malfunctioning.
- Only use the attachments on the recommended muscle groups.
- Do not use on injured or swollen areas, or if you have open wounds, bruises, broken bones, or skin conditions.
- Do not use for more than 60 seconds on a single area. Move the device slowly across muscles.
- If blood pressure is a concern, please check your blood pressure before and after percussive therapy, as vibrations increase blood circulation.
- If you suffer from a heart condition, please consult your medical provider before engaging in percussive therapy.

- If you are pregnant, please consult your medical provider before engaging in percussive therapy.

Safety

- Stop use immediately if you feel pain, discomfort, dizziness, or notice any abnormal reaction.
- Do not disassemble or attempt to repair the device yourself. Contact customer service for assistance.
- Do not operate while sleeping, drowsy, or under the influence of drugs or alcohol.
- The product should be used by adults only and kept out of reach of children.
- Consult your doctor before use if you have any medical conditions, such as heart problems, diabetes, or if you are pregnant.
- Do not use the massage gun while it is charging.
- Do not put objects or other things in the Massage Head opening hole of the massage gun.

Troubleshooting

Question	Causes	Solutions
Product stopped working	Dead battery	Recharge the device
	Pressure applied is too high	Ease off pressure
	Device is on charge	Unplug device
Lasts only few minutes even after full charge	Dead battery	Contact service centre for support
	Insufficient charging time	Charge the unit for at least 4 hours
	Charging cord damaged	Try charging the device with new charging cable
	Adapter broken	Try a new adapter
Massage heads popping out	Incorrect insertion	Push massage head all the way in before use
Warm head isn't hot	A speed level hasn't been selected	Select a speed level and wait for a few seconds

Environmental Protection and Recycling

Safe Appliance Disposal:

At the end of its working life, the battery must not be disposed of as urban waste. It must be taken to a special authority differentiated waste collection centre or to a dealer providing

this service. Disposing of a household appliance separately avoids possible negative consequences for the environment and health deriving from inappropriate disposal. This enables the constituent materials to be recovered to obtain significant savings in energy and resources.

As mentioned, the built-in rechargeable battery in this product contains substances that could be harmful to the environment. To ensure the responsible disposal and recycling of the device, please follow your country's rules for the separate collection of electrical products and batteries, or refer to the following steps:

1. Activate the device's power switch and let it run until it stops functioning, ensuring the battery is entirely drained.
2. Before disassembling the device and removing its battery, always ensure the device is powered off and disconnected from any power sources.
3. Once the battery has been extracted, never reconnect the device to its charger. Additionally, be cautious to avoid causing a short circuit between the battery's positive and negative terminals.
4. Then follow local authority guidelines to dispose of the removed battery.

Customer Service

For assistance, please visit our website: [OOZZ.com](https://www.oozz.com). You can also call or WhatsApp us on 0800 193 6789 (UK), +44800 193 6789 (International) or email us at care@oozz.com.

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