oozz

0Z-3621


Importer: Elite Capital Limited

Address: Studio 4 Clarks Courtyard, 145 Granville Street,
Birmingham, B1 1SB, United Kingdom.

Website: oozz.com

Email: care@oozz.com

OZ-3621

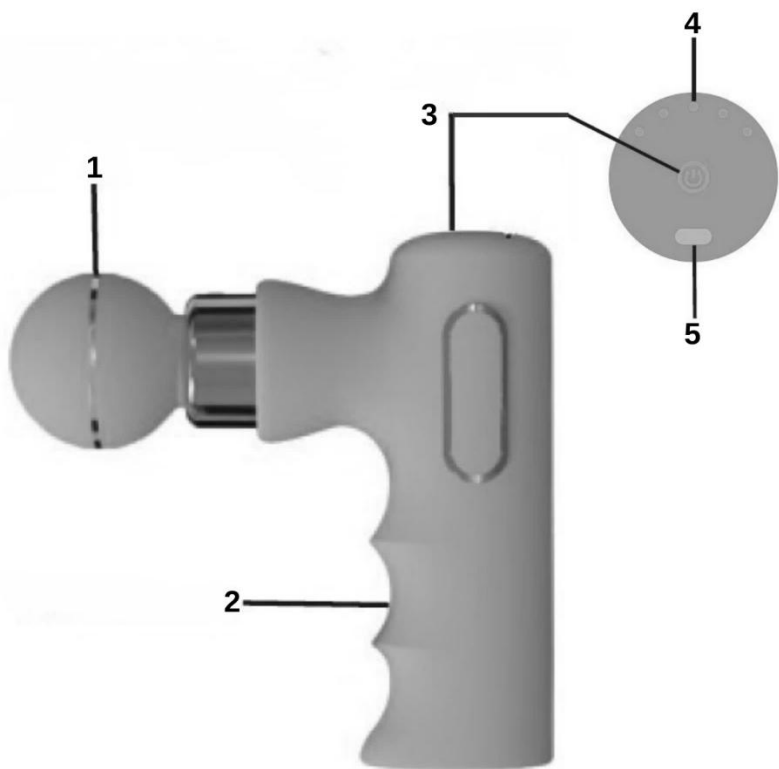
6

7

8

9

10

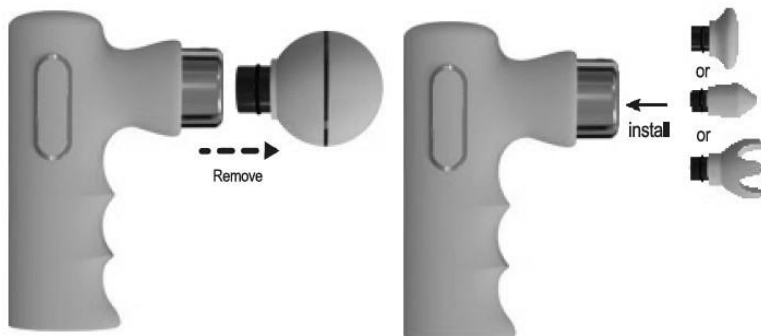Welcome to OZZZ and congratulations on your new purchase! Secure your peace of mind by activating your 12-month warranty and claiming extra 6 months for free at www.oazz.com/warranty. Beyond our warranty, you are also covered by your country's legal consumer rights, ensuring complete protection for your purchase.

1. Massage Head
2. Handle
3. Power Button
4. Indicator Lights
5. Charging Interface
6. Charging Cable
7. Flat Head (Suitable for all areas of the body)
8. Sphere Head (Ideal for large muscle groups)
9. Cylinder Head (For joints and targeted deep-tissue areas)
10. U-Head (For the spine and the area behind the ankle)

Product Details

Product Name	OZ-3621 Mini Massage Gun
Model Number	OZ-3621
Charging Input	DC 5V 1A
Charging Time	≤4 Hours
Usage Time	80 mins (Max level)

Charging Instructions:

- Make sure the device is switched OFF before charging.
- Charge the product fully before you use it for the first time and after a long period of disuse. For optimal performance, it is best to use the device when fully charged.
- When the power is low, the indicator light will flash to remind you to charge the device.
- The charging indicator will blink while charging and turn solid when fully charged.

How to use:

- Long press the power button to Turn ON or Turn OFF the device.
- The power button has 5 intensity levels. It starts at level 1, and each press increases the level by one.

Precautions:

- Charge only with DC 5V 1A certified adapter.
- Do not overcharge or charge near water. Do not expose to water or moisture. This device is not waterproof.
- Do not use if damaged or malfunctioning.
- The product should be used by adults only and keep out of reach of children.
- For your safety, do not use the massage gun while it is charging.

- Use the device only on soft tissues of the body without causing pain or discomfort. Do not use it on head, bones of sensitive parts of the body.
- Stop using the device once there is any sign of pain or discomfort.
- Only use the attachments on the recommended muscle groups.
- To protect the bones and muscles, avoid long-term fixed-point massage.
- Do not immerse the device in water or allow water to enter the device.
- Do not use on injured or swollen areas, or if you have open wounds, bruises, or skin conditions.
- Consult your doctor before use if you have any medical conditions, such as heart problems, diabetes, or if you are pregnant.
- Do not use for more than 60 seconds on a single area. Move the device slowly across muscles.
- Do not operate while sleeping, drowsy, or under the influence of drugs or alcohol.
- Use only the provided charger and accessories.
- Stop use immediately if you feel pain, discomfort, dizziness, or notice any abnormal reaction.
- Do not disassemble or attempt to repair the device yourself. Contact customer service for assistance.

Cleaning:

- Power off and disconnect before cleaning.
- To clean, wipe the surface of the device with a damp towel and dry it with a soft towel.

Customer Service

For assistance, please visit our website: OOZZ.com. You can also call or WhatsApp us on 0800 193 6789 (UK), +44800 193 6789 (International) or email us at care@oozz.com.

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