oozz

0Z-3631




Importer: Elite Capital Limited

Address: Studio 4 Clarks Courtyard, 145 Granville Street,
Birmingham, B1 1SB, United Kingdom.

Website: oozz.com

Email: care@oozz.com

OZ-3631

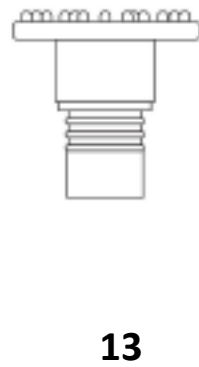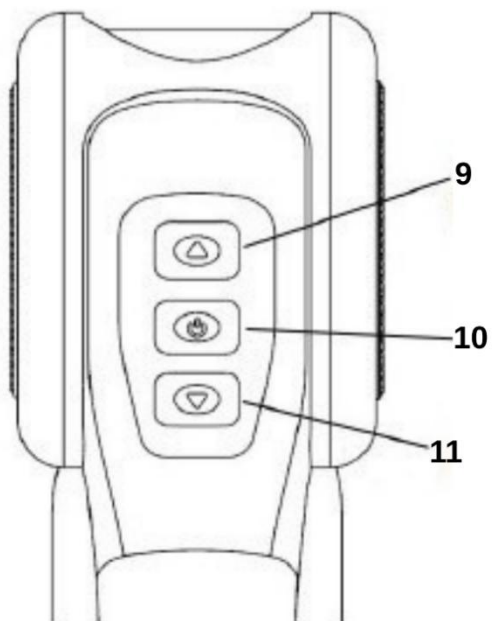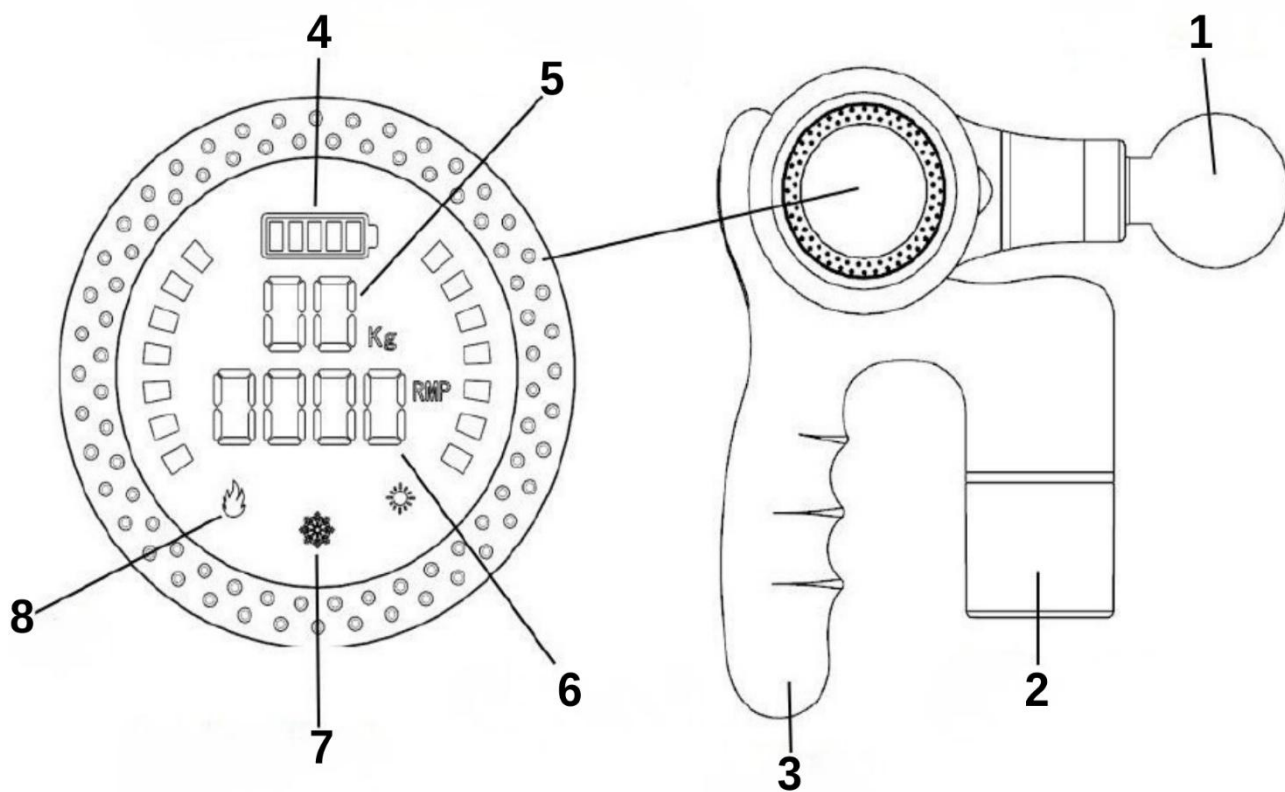14

15

16


17

18


19

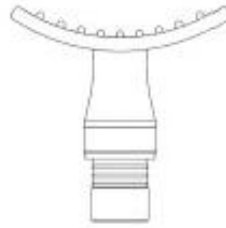20

21

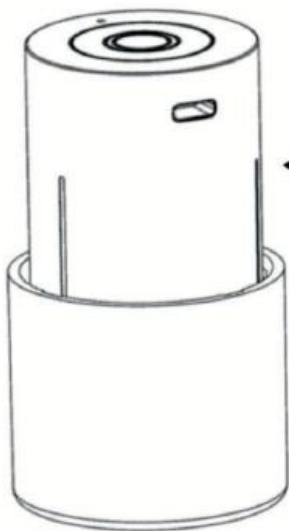Type C

USB

(Self Provided)
Note: The Output
is 5V-1A mobile
phone charging head

Welcome to OOOZ and congratulations on your new purchase! Secure your peace of mind by activating your 12-month warranty and claiming extra 6 months for free at www.oozz.com/warranty. Beyond our warranty, you are also covered by your country's legal consumer rights, ensuring complete protection for your purchase.

1. Massage Head
2. Battery Pack
3. Handle
4. Battery Indicator
5. Force Indicator
6. Speed Indicator
7. Cold Compress Indicator
8. Hot Compress Indicator
9. Speed Increase Button
10. Power Button
11. Speed Decrease Button
12. Charging Cable
13. Flat Head (Suitable for large muscle groups such as the back, waist, arm)
14. Sphere Head (For large muscle groups across the body)
15. Cylinder Head (For joints and targeted deep-tissue areas)
16. U-Head (For the spine and area behind the ankle)
17. Finger Head (For deep-tissue area such as soles, joints, palms)
18. Silicone Swamm Head (for the whole body, especially medium to large muscle areas, due to its soft and gentle surface)
19. Round Scoop Head (For small muscle areas)
20. Camber Head (Ideal for the legs, arms, and neck)
21. Hot & Cold Head (Features a heating function that reaches up to 50°C to relax muscles, and a cooling function that drops to 10-15°C to reduce soreness and swelling)

Product Details

Product Name	OZ-3631 Massage Gun
Model Number	OZ-3631
Charging Input	DC 5V 2A
Speed (Maximum and Minimum)	1600-3000 r/min
Charging Time	≤5 Hours
Usage Time	120 mins

Charging Instructions:

- Make sure the device is switched OFF before charging.
- Charge the product fully before you use it for the first time and after a long period of disuse. For optimal performance, it is best to use the device when fully charged.

- To charge the device, remove the battery pack from the massage gun. Connect the Type-C converter to the battery, then plug in the USB cable and connect it to a power source or USB adapter.
- Once the battery is fully charged, attach the battery back onto the device.
- When the device power is low, the Battery Indicator will flash to remind you to charge the device.
- When the battery pack is charging, the indicator light on the battery pack will turn red. It will change to Green when fully charged.

Remark: This massage gun includes an extra replaceable battery pack, so you won't run out of power during your massage sessions.

How to use:

- Long press the power button to Turn ON or Turn OFF the device.
- Press '▲and▼' button to adjust massage force (1600-3000). The LED Display will correspond to the selected level of force.
- Click the '▲' button to start the vibration and increase the force intensity and click '▼' key to decrease the force intensity.
- **Hot & Cold Head:** Press the start button once to activate the cold (blue) mode and press it again to switch to the hot (red) mode.

Note: To ensure optimal performance, set the massage gun to **the lowest speed** when using the hot & cold head.

Precautions:

- Charge only with DC 5V 2A certified adapter.
- Do not overcharge or charge near water. Do not expose to water or moisture. This device is not waterproof.
- Do not use if damaged or malfunctioning.
- The product should be used by adults only and keep out of reach of children.
- For your safety, do not use the massage gun while it is charging.
- Use the device only on soft tissues of the body without causing pain or discomfort. Do not use it on head, bones of sensitive parts of the body.
- Stop using the device once there is any sign of pain or discomfort.
- Only the attachments on the recommended muscle groups.
- To protect the bones and muscles, avoid long-term fixed-point massage.
- Do not immerse the device in water or allow water to enter the device.
- Do not use on injured or swollen areas, or if you have open wounds, bruises, or skin conditions.
- Consult your doctor before use if you have any medical conditions, such as heart problems, diabetes, or if you are pregnant.
- Do not use for more than 60 seconds on a single area. Move the device slowly across muscles.
- Do not operate while sleeping, drowsy, or under the influence of drugs or alcohol.
- Use only the provided charger and accessories.
- Stop use immediately if you feel pain, discomfort, dizziness, or notice any abnormal reaction.

- Do not disassemble or attempt to repair the device yourself. Contact customer service for assistance.

Cleaning:

- Power off and disconnect before cleaning.
- To clean, wipe the surface of the device with a damp towel and dry it with a soft towel.

Customer Service

For assistance, please visit our website: OOZZ.com. You can also call or WhatsApp us on 0800 193 6789 (UK), +44800 193 6789 (International) or email us at care@oozz.com.

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