# **oozz**

## **OZ-7213**


Importer: Elite Capital Limited

Address: Studio 4 Clarks Courtyard, 145 Granville Street,  
Birmingham, B1 1SB, United Kingdom.

Website: [oozz.com](http://oozz.com)

Email: [care@oozz.com](mailto:care@oozz.com)


OZ-7213

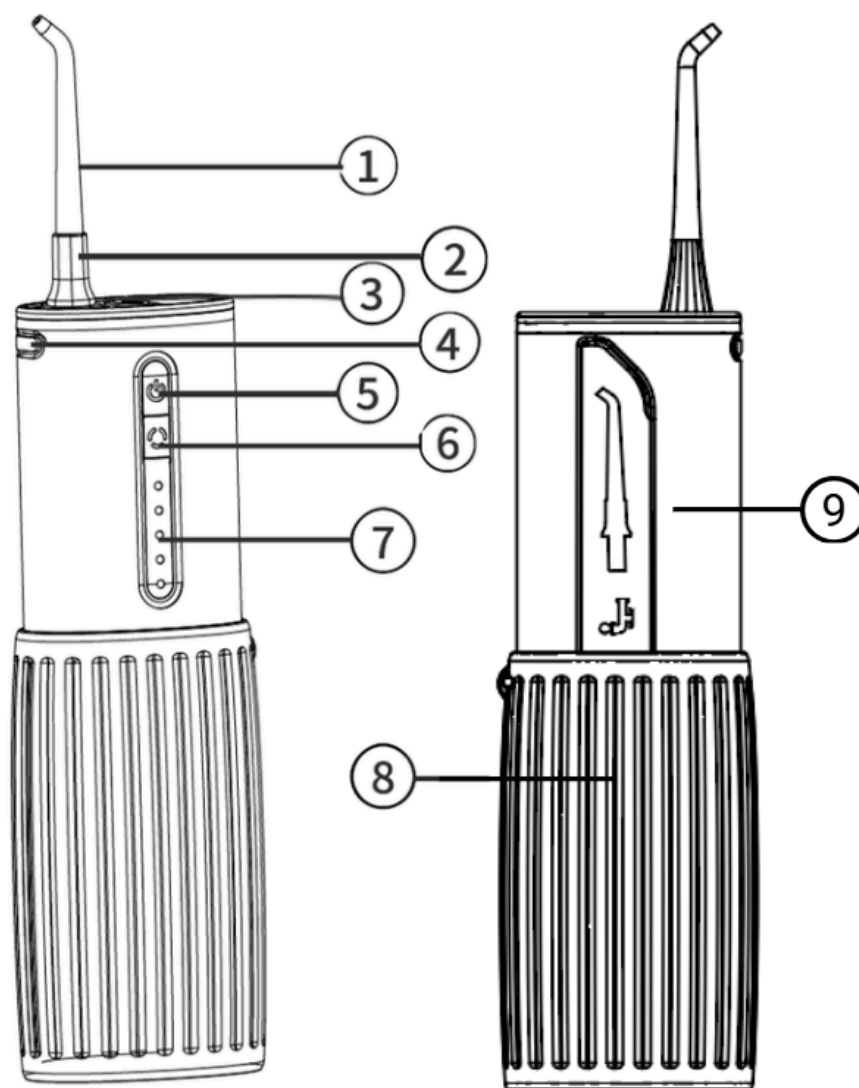
10

11

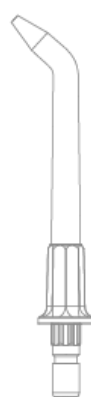12

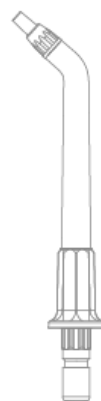13

14

Welcome to OZZZ and congratulations on your new purchase! Secure your peace of mind by activating your 12-month warranty and claiming an extra 6 months for free at <https://www.oazz.com/pages/warranty>. Beyond our warranty, you are also covered by your country's legal consumer rights, ensuring complete protection for your purchase.

- 1. Nozzle
- 2. Nozzle Rotation Dial
- 3. USB Charging Port
- 4. Nozzle Release Button
- 5. Power Button
- 6. Mode Selection Button
- 7. LED Indicator
- 8. Water Tank
- 9. Nozzle Storage Room / Water Injection
- 10. USB Charging Cable
- 11. Standard Tip X3
- 12. Orthodontic Tip X1
- 13. Periodontal Tip X1
- 14. Tongue Scraper Tip X1

Product Details

|                     |                              |
|---------------------|------------------------------|
| Product Name        | OZ-7213 Travel Water Flosser |
| Model Number        | OZ-7213                      |
| Charging Input      | 5V 1.0A                      |
| Water Tank Capacity | 150ML                        |
| Charging Time       | 5 Hours                      |
| Usage Time          | Up to 120 days               |

Nozzle Tips: Function and Usage

| Nozzle Type        | Usage                                                                         | Recommended Replacement Interval |
|--------------------|-------------------------------------------------------------------------------|----------------------------------|
| Standard Tip       | Ideal for everyday cleaning and removing food debris between teeth and gums.  | Every 3 Months                   |
| Orthodontic Tip    | Designed for cleaning around braces, brackets, and other orthodontic devices. |                                  |
| Periodontal Tip    | Gently cleans deep below the gumline to treat gum disease and improve health. |                                  |
| Tongue Scraper Tip | Helps remove bacteria and debris from the tongue, promoting fresher breath.   |                                  |

## Instructions For Use

### 1. On/Off Button:

- **Start Flossing:** Insert the tip, then press the power button to begin irrigating. The device has a 2-minute auto-stop timer or can be turned off manually by pressing the power button again.
- **Turn Off:** Once finished, press the power button to turn off the device.
- **Empty and Drain:** Empty any remaining water from the tank, then run the device for 2-3 seconds to drain residual water.

### 2. Inserting and Removing the Nozzle:

- Insert the nozzle into the top of the water flosser handle, pressing firmly until you hear a click, confirming it is securely in place (Figure 1).

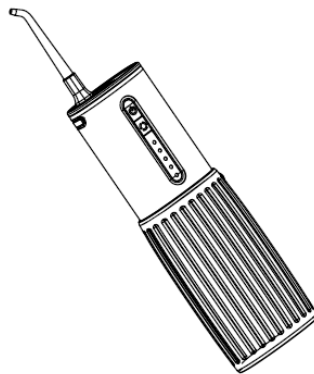
Figure 1

- To remove the nozzle, press the tip release button and gently pull the nozzle out (Figure 2).

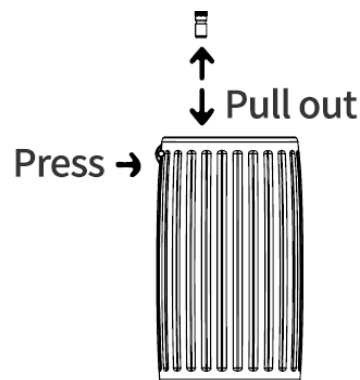
Figure 2

### Note:

- Do not press the tip release button while the device is operating.
- If the nozzle head doesn't spray water immediately after turning on the device (particularly after first use or long periods of inactivity), follow these steps:
  1. Add a small amount of water to the tip insertion area.
  2. Insert a new tip and press the power button to activate the device.

### 3. Installing and Filling the Water Tank:

- Open the water tank cap, fill with water, and close the cap securely (Figure 3).

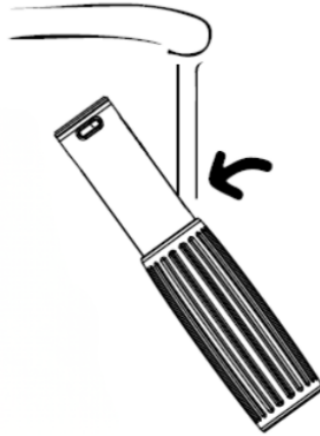
Figure 3

- To detach or attach the water tank, simply twist it.
  - Place the device on a flat surface.
  - Hold the water tank with your left hand and twist the upper part to the right with your right hand to detach.
  - To reattach, twist in the opposite direction.

**Note:** Use water or a dental professional-recommended solution only when filling the tank.

### 4. Pressure Mode Settings:

This water flosser offers 5 pressure modes, easily adjustable by pressing the “” Mode button. The pressure settings are as follows:

- **Gentle:** For beginners
- **Medium:** Daily cleaning
- **Strong:** Removes stubborn stains
- **Pulse:** Alternating strong/weak flow
- **Custom:** Long press the Mode button to adjust water flow

**Note:** Check the set mode with the corresponding mode indicator light.

### 5. How to Charge:

Connect the provided USB cable directly to the device. Then plug the USB into a 5V1A USB adapter. The product must be charged via a laptop/computer/power bank USB port (Figure 4). After charging remove the power plug from the USB outlet and then remove the USB out of the product.

- Charge the product fully before you use it for the first time and after a long period of disuse.
- Make sure the device is switched OFF before charging.

- The LED light will light up, confirming the charging process has begun.
- A full charge takes 5 hours providing up to 120 days of operation. For optimal performance, it's best to use the device when fully charged.

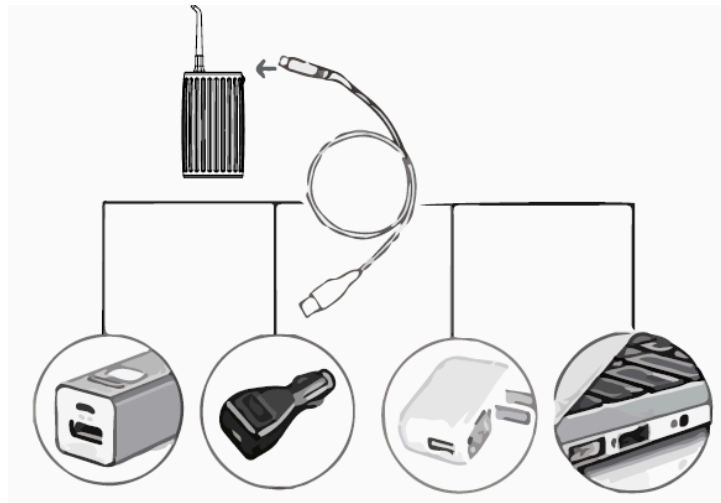
Figure 4

### Cleaning & Maintenance

Please note: Although this appliance is waterproof (IPX8 Rating), do not immerse the device in water or any other liquid. To maintain optimal performance and longevity, make sure that the product is cleaned and thoroughly dried before its next use.

- Before cleaning, ensure the device is switched OFF and disconnected from any charging cables.
- Empty the water tank after each use and run the device for 2-3 seconds to drain any remaining water to prevent bacteria buildup.
- Always use only water or neutral cleaning agents to clean the appliance. Never use caustic or abrasive agents (e.g. vinegar or lime scale removers), since they may cause damage to the appliance.
- Keep the appliance away from high temperature and avoid direct sunlight.
- Do not submerge the device in water or any other liquid.
- Wipe the device with a soft cloth and a mild, non-abrasive cleaner when necessary.
- Clean the water tank with soapy water; use a brush to remove any dirt inside and wipe clean with a cloth.
- Ensure the water tank is completely dry before storing the device for long periods.

### Precautions

#### Charging:

- Ensure the device is powered off before connecting the charging cable. Avoid forcefully pulling on the power supply wire or the device.
- Always attach the power cord to the product first, then to the power outlet.

- If the USB cable accidentally falls into water, unplug the adapter first before retrieving the product.
- Do not unplug the adapter with wet hands.
- Do not plug this device into a voltage system that is different from the voltage system specific on the device or charger.
- Prevent overcharging by unplugging the device once the designated charging time is reached.
- Always keep the product away from water sources while charging, and avoid charging it near baths, sinks, or any liquid containers.
- Do not charge on soft surfaces like beds and bedding, especially in enclosed spaces. Should the product overheat during charging, disconnect it immediately.
- Avoid charging in high-temperature areas or near open flames to mitigate fire risks.
- Avoid damaging, breaking, or modifying the USB cable. Do not stretch, bend or twist the USB cable excessively, and do not place heavy objects or trap it between objects.
- If the power cord becomes damaged, replace it using only compatible cords or components purchased directly from the original manufacturer.

#### **Usage:**

- This product is designed solely for cleaning your teeth, gums and tongue. Do not use it for any other purpose.
- Do not use the product if you lose the nozzle.
- Only use tips and accessories that are recommended.
- Do not press the nozzle too firmly against your teeth or gums, or it may cause injury.
- People suffering from tooth or gum pain may not be able to use the appliance and should consult their dentist before use.
- Do not insert any foreign object into any opening on the product.
- Do not use the product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and shorten lifespan of the product. Use of these products may void warranty.
- Remove any oral jewellery prior to use of this product.
- Always power off the device prior to cleaning or performing maintenance.
- If the device experiences any malfunction or damage, immediately power it off and avoid unauthorised disassembly.
- This water flosser is waterproof (IPX8 Rating) however do not immerse the product in water to prevent malfunctioning, fire, electric shock or explosion.
- Ensure the product is completely dry before storing to maintain its condition and functionality.
- We recommend avoiding storage in damp, humid areas and high temperatures. Do not throw product into the fire or heat it.

**Safety:**

- The appliance can be used by children above 8 years of age however supervision is recommended to ensure safe operation. Always store the device out of reach of unsupervised children to prevent any accidental harm.
- Children should not play with the device as this appliance is not a toy.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. Also, the deadly amoeba, naegleria fowleri, may be present in some tap water or non chlorinated well water and may be fatal if directed in these areas.
- Keep the appliance, cable, nozzle out of reach of children under 8 years. The nozzle is small enough to be swallowed and can be lodged in the throat.

**Troubleshooting**

| Question                                      | Causes                                                                                         | Solutions                                                                                           |
|-----------------------------------------------|------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------|
| Water tank leaks                              | Seal ring is abnormal                                                                          | Remove the seal ring and reinstall                                                                  |
| Low water pressure                            | Low battery                                                                                    | Recharge the device                                                                                 |
|                                               | Nozzle head blocked by dirt or debris                                                          | Replace with new nozzle head                                                                        |
| Lasts only few minutes even after full charge | Dead battery                                                                                   | Contact service centre for support                                                                  |
|                                               | Insufficient charging time                                                                     | Charge the unit for at least 5 hours                                                                |
|                                               | Charging cord damaged                                                                          | Try charging the device with new charging cable                                                     |
|                                               | Adapter broken                                                                                 | Try a new adapter                                                                                   |
| Motor makes sound but water does not come out | Water tank is empty                                                                            | Fill the water tank                                                                                 |
|                                               | Suction valve dislocation                                                                      | Keep device on, fill water to the hole where nozzle head goes into                                  |
|                                               | Nozzle head blocks                                                                             | Replace with new nozzle head                                                                        |
| Motor doesn't make sound                      | Low battery                                                                                    | Activate the battery by charging the unit for 16 hours                                              |
| Product stopped working suddenly              | Overheating protection is activated (continuously used 5 tanks of water with highest pressure) | Excessive heat will cause the machine to stop. Switch off the device, cool for 2 hours, and restart |

## Environmental Protection and Recycling

### Safe Appliance Disposal

At the end of its working life, the battery must not be disposed of as urban waste. It must be taken to a special authority differentiated waste collection centre or to a dealer providing this service. Disposing of a household appliance separately avoids possible negative consequences for the environment and health deriving from inappropriate disposal. This enables the constituent materials to be recovered to obtain significant savings in energy and resources.

As mentioned, the built-in rechargeable battery in this product contains substances that could be harmful to the environment. To ensure the responsible disposal and recycling of the device, please follow your country's rules for the separate collection of electrical products and batteries, or refer to the following steps:

1. Activate the device's power switch and let it run until it stops functioning, ensuring the battery is entirely drained.
2. Before disassembling the device and removing its battery, always ensure the device is powered off and disconnected from any power sources.
3. Once the battery has been extracted, never reconnect the device to its charger. Additionally, be cautious to avoid causing a short circuit between the battery's positive and negative terminals.
4. Then follow local authority guidelines to dispose of the removed battery.

### Customer Service

For assistance, please visit our website: [OOZZ.com](http://OOZZ.com). You can also call or WhatsApp us on 0800 193 6789 (UK), +44800 193 6789 (International) or email us at [care@oozz.com](mailto:care@oozz.com).

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