oozz

0Z-7213


Importer: Elite Capital Limited

Address: Studio 4 Clarks Courtyard, 145 Granville Street,
Birmingham, B1 1SB, United Kingdom.

Website: oozz.com

Email: care@oozz.com

OZ-7213

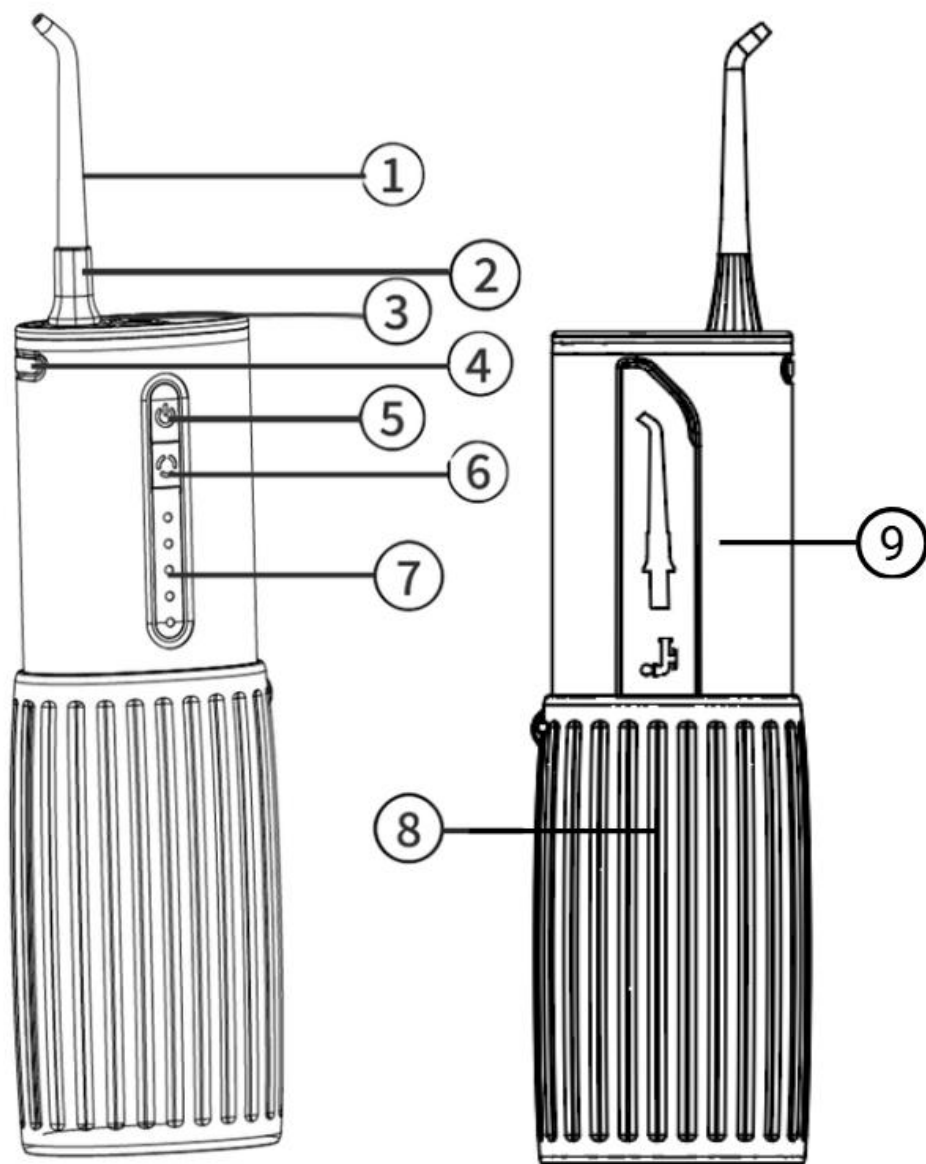
10

11






12

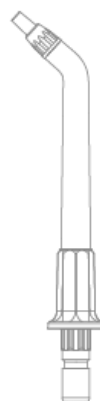13

14

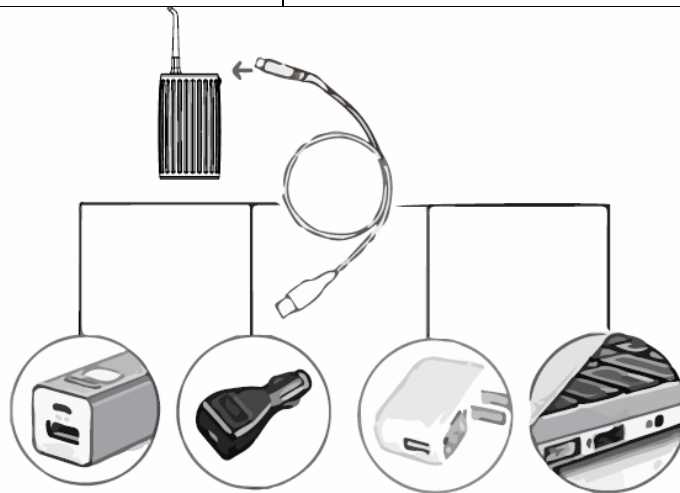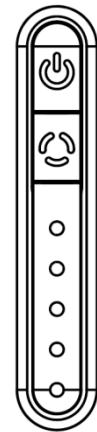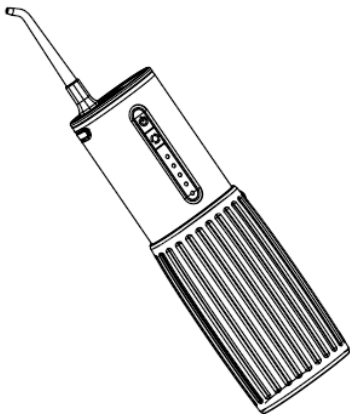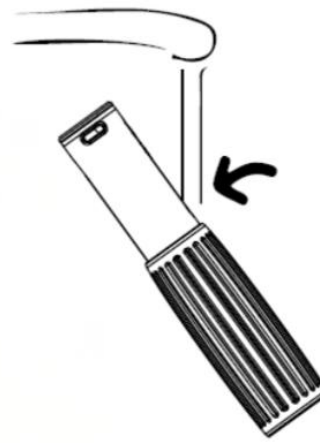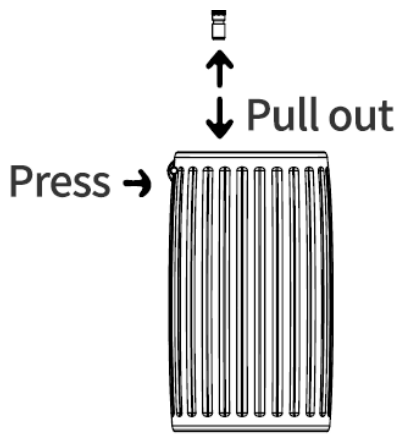Welcome to OOZZ and congratulations on your new purchase! Secure your peace of mind by activating your 12-month warranty and claiming an extra 6 months for free at www.oozz.com/warranty. Beyond our warranty, you are also covered by your country's legal consumer rights, ensuring complete protection for your purchase.

1. Nozzle
2. Nozzle Rotation Dial
3. USB Charging Port
4. Nozzle Release Button
5. Power Button
6. Mode Selection Button
7. LED Indicator
8. Water Tank
9. Nozzle Storage Room / Water Injection
10. USB Charging Cable
11. Standard Tip X3
12. Orthodontic Tip X1
13. Periodontal Tip X1
14. Tongue Scraper Tip X1

Product Details

Product Name	OZ-7213 Travel Water Flosser
Model Number	OZ-7213
Charging Input	5V 1.0A
Water Tank Capacity	150ML
Charging Time	5 Hours
Usage Time	Up to 120 days

Nozzle Tips: Function and Usage

Nozzle Type	Usage	Recommended Replacement Interval
Standard Tip	Ideal for everyday cleaning and removing food debris between teeth and gums.	Every 3 Months
Orthodontic Tip	Designed for cleaning around braces, brackets, and other orthodontic devices.	
Periodontal Tip	Gently cleans deep below the gumline to treat gum disease and improve health.	
Tongue Scraper Tip	Helps remove bacteria and debris from the tongue, promoting fresher breath.	

Charging Instructions:

- Make sure the device is switched OFF before charging.
- Charge the product fully before you use it for the first time and after a long period of disuse. For optimal performance, it is best to use the device when fully charged.
- While charging, 5 LED lights flash from bottom to top. All 5 lights stay on when fully charged.
- Battery Level: All 5 lights mean full power, 3 lights indicate half, 1 light shows low power, and a flashing light means it's time to charge.
- Connect provided USB cable to device and plug into a 5V 1A USB adapter or power source (laptop/computer/power bank).

How to Use:

- Insert the tip, then press the power button to begin irrigating; once finished, press the power button to turn off the device.
- The device has a 2-minute auto-stop timer or can be turned off manually by pressing the power button again.
- Pressure Mode Settings. This water flosser has 5 modes:
 - Gentle – for beginners
 - Medium – for daily cleaning
 - Strong – removes stubborn stains
 - Pulse – alternating strong/weak flow
 - Custom – long press the Mode button to adjust water flow
- To remove the nozzle, press the nozzle tip release button and gently pull the nozzle out.

Precautions:

- Charge only with 5V 1A certified adapter.
- Do not overcharge or charge near water.
- Do not use if damaged or malfunctioning.
- Store between 5°C and 35°C. Avoid damp areas.
- Keep out of reach of unsupervised children under 8 years of age.

Cleaning:

- Power off and disconnect before cleaning.
- Empty the water tank after each use and run the device for 2-3 seconds to drain any remaining water and prevent bacteria buildup.
- Ensure the water tank is completely dry before storing the device for long periods.

Customer Service

For assistance, please visit our website: OOZZ.com. You can also call or WhatsApp us on 0800 193 6789 (UK), +44800 193 6789 (International) or email us at care@oozz.com.

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